




Apollo CRO

150° Color, AF and FFA
Video Recording
Offline Retina Ass<u>istant</u>

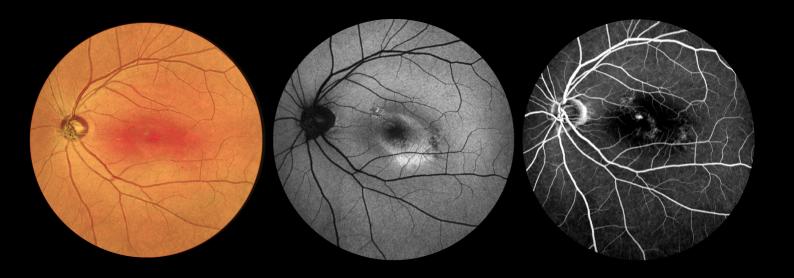

Microclear Apollo ultra-widefield CRO with AF and FFA

No compromises Fundus Fluoroangiography videos and Autofluorescence



Autofluorescence

Single shot or combine multiple shots to increase the resolution



Fluorescein Angiography

Maximum control of your Angiography examination with all the features you need



Video Recording

Don't miss the initial stages and store and review from any workstation



Pupil size 1.5mm

Optimal photos even on really small pupils resulting on higher success rate of gradable images



150° one shot and up to 240° in two shots

Capture two to six photos and combine in a panorama view easily thanks to the fast and precise joystick control



Optical Zoom 150° - 90° - 45° - 22° without changing lens

One click switch Field of View to zoom into areas of interest



Strong penetration: cataract free

Strong Confocal design allows to separate better the retina layers and to penetrate stronger through cataract and eye opacities.



Unlimited offline Retina Assistant included

Ready to use Risk Assessment software with 19 pathologies detection to support your decisions




Cells and blood flow tools option


Optional module available for 5° cellular view imaging with software for retinal cell count, blood flow velocity and vascular thickness measurement

Full fledged Autofluorescence and Fluorescein Angiography

Capture the big picture with 150° photos, and one-touch zoom to 90°, 45° and 22° views. Don't miss FFA early stages with unlimited videos at any angle combination.

Apollo CRO is fully configurable and has all the flexibility required to empower your FFA examination, yet easy to use with a fast learning curve.

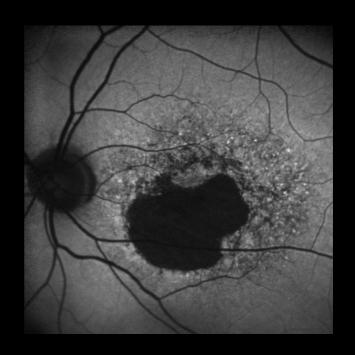

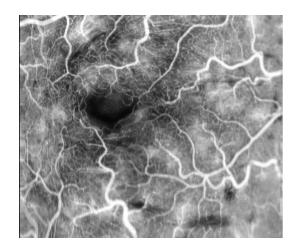

Upgrade your imaging capabilities

Obtain ultra-widefield images with easy and fast capture of 150° retinal photos in all modalities: Infrared, Redfree, Colour, Autofluofluorescence, Fundus Fluoroangiography, and create panorama views up to 240°.

Controlling Apollo feels comfortable thanks to the fast and precise joystick that allows to take a good photo in a matter of a few seconds.

The Strong Confocal optical design delivers a high success rate of gradable images even with small pupils and adverse light situations.







Apollo CRO technical specifications

Imaging Modalities	Infrared - 785 nm Redfree - 488 nm Color Autofluorescence Fluorescein Angiography photo and video
Field of View (eye center)	One shot single capture Optical zoom 150° / 90° / 45° / 22°
Automatic Montage	Up to 240° mosaic 2 to 9 images
Non Mydriatic	Minimum pupil 1.5mm
Resolution	5 microns
Fixation Targets	Internal 1 central and 8 peripheral fixation targets External flexible target
Diopter Compensation	-15D - +15D
Software Included	Patient management Retina review and annotation Multiple review stations DICOM compatible
Retina Assistant	19 pathologies detection Offline, no internet requried Unlimited license use
Hardware Included	PC computer Monitor Table (printer)
Power	Voltage: 100-240V~50-60Hz

Microclear is committed to the highest standard of retina imaging

Suzhou MicroClear Medical Instruments Co., Ltd

Tel: +86-512-67067163 Fax: +86-512-67066173

E-mail: marketing@microcleartech.com

Suite 1601-1602, G2 Building, No.88 Jinjihu Avenue, Suzhou Industrial Park, Suzhou, Jiangsu P.R.China